# VenaTech® Convertible™

## Vena Cava Filter System

es Sistema de filtro de Vena Cava

en Instructions for Use Jugular Approach

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en Manufactured for:

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**B. Braun Interventional Systems Inc.**

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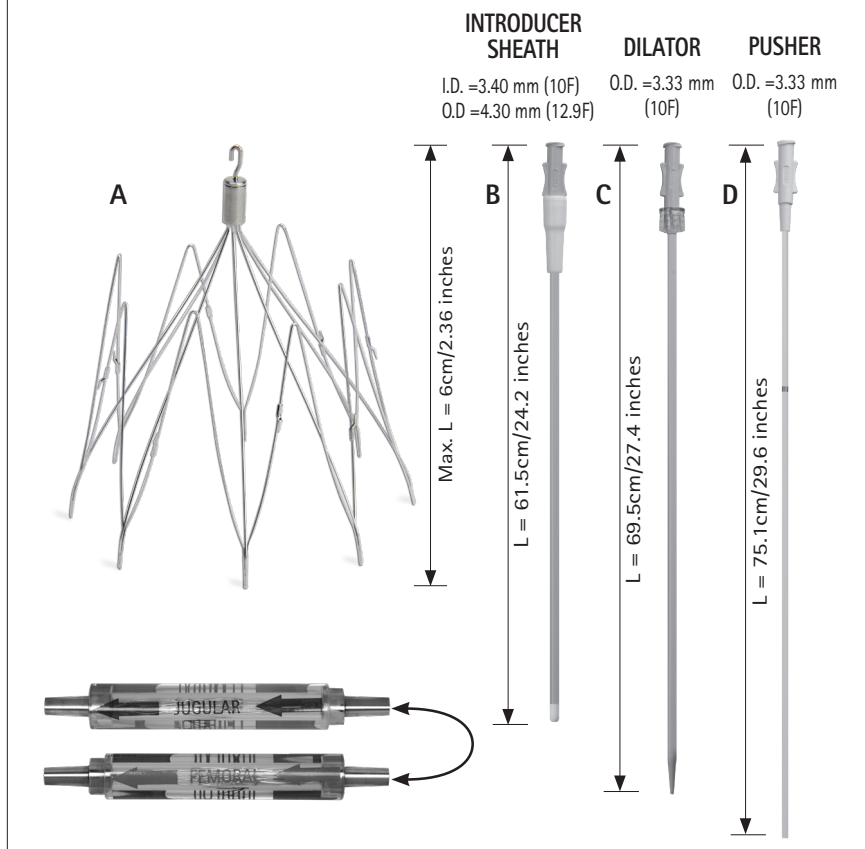
**Fig. A** en 1 VenaTech® Convertible™ Vena Cava Filter pre-loaded  
 in its cartridge for use with femoral or jugular approaches.

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Indications ..... 5

Contraindications..... 6

Potential Adverse Events / Technical Complications..... 6

Important Information, General Warnings and Precautions ..... 7-10

MRI Safety Information..... 10

Clinical Study Summary ..... 11-17

Implantation Procedure..... 18-21

Conversion Procedure ..... 22-25

Final Documentation ..... 26

Conversion with Accessories..... 26

**WARRANTY:** B. BRAUN INTERVENTIONAL SYSTEMS INC. WARRANTS THAT THIS PRODUCT WAS MANUFACTURED ACCORDING TO APPLICABLE STANDARDS AND SPECIFICATIONS. PATIENT CONDITION AND CLINICAL TREATMENT MAY AFFECT THE PERFORMANCE OF THIS PRODUCT. USE OF THIS PRODUCT SHOULD BE IN ACCORDANCE WITH THE INSTRUCTIONS PROVIDED AND AS DIRECTED BY THE PRESCRIBING PHYSICIAN.

## **DEVICE DESCRIPTION**

**CAUTION: Federal law (USA) restricts this device to sale by or on the order of a physician.**

The VenaTech® Convertible™ Vena Cava Filter System consists of an implantable vena cava filter and the introducer accessories. The VenaTech Convertible Vena Cava Filter is a flexible, symmetrical, cobalt-chromium metal alloy, self-expanding vena cava filter intended for permanent implantation into the infrarenal inferior vena cava to provide protection against pulmonary embolism. When clinically indicated, the clot-trapping features of the VenaTech Convertible Vena Cava Filter can be deactivated by percutaneously converting the filter to an open configuration.

The VenaTech Convertible Vena Cava Filter is pre-loaded in a cartridge and provided as a system with introducer accessories; system components include:

- VenaTech Convertible Vena Cava Filter (pre-loaded in a cartridge)
- Introducer System:
  - Introducer Sheath
  - Dilator
  - Pusher

The VenaTech Convertible Vena Cava Filter System is supplied STERILE and non-pyrogenic in a peel-open package. The system is sterilized using Ethylene Oxide. The system is intended for single-use only. Dispose of system components after filter implantation. Components may not be cleaned, re-sterilized, or reused. Store at room temperature. Avoid freezing and excessive heat. Do not remove from sterile package until ready for use.

## **OTHER REQUIRED ACCESSORIES**

Accessories required for the filter implantation procedure pre-placement cavogram but not provided in the package:

- Seldinger needle (when Seldinger technique is used)
- 'J' Guidewire (0.035 in / 0.89 mm, 70 cm length)
- Catheter introducer sheath (9F I.D.)
- Catheter (angiography catheter or other catheter of choice)
- Scalpel

**Accessories required for the filter implantation procedure but not provided in the package:**

- Coated 'J' Guidewire (0.035 in / 0.89 mm, 150 cm length)

**Accessories required to convert the filter to an open configuration but not provided in the package:**

Physician must have the following items available to facilitate conversion and ensure filter legs are fully apposed to vessel wall:

- 'J' Guidewire (0.035 in / 0.89 mm, 150 cm length)
- Catheter introducer sheath (9F I.D.)
- Amplatzer GooseNeck™ snare (size matched to the diameter of the vena cava) and catheter
  - 15 mm loop and 6F catheter (I.D. 1.6 mm), length ≥ 102 cm
  - 20 mm loop and 6F catheter (I.D. 1.6 mm), length ≥ 102 cm
  - 25 mm loop and 6F catheter (I.D. 1.6 mm), length ≥ 102 cm
- Accessory device such as a diagnostic catheter, guidewire, or angioplasty balloon

**INDICATIONS**

The VenaTech® Convertible™ Vena Cava Filter System is indicated for the prevention of recurrent pulmonary embolism via placement in the vena cava in the following situations:

- Pulmonary thromboembolism when anticoagulants are contraindicated;
- Failure of anticoagulant therapy in thromboembolic diseases;
- Emergency treatment following massive pulmonary embolism where anticipated benefits of conventional therapy are reduced; and
- Chronic, recurrent pulmonary embolism where anticoagulant therapy has failed or is contraindicated

When clinically indicated after implantation, The VenaTech Convertible Vena Cava Filter can be converted to an open configuration to discontinue filtration according to the 'Conversion Procedure' in the Instructions for Use.

GooseNeck is a trademark of Covidien.

## **CONTRAINDICATIONS**

The VenaTech® Convertible™ Vena Cava Filter System is contraindicated for the following situations:

- Vena cava filters must not be implanted in patients with a risk of septic embolism.
- Patients whose vena cava diameter is greater than (>) 28 mm (due to the risk of device migration).
- Pregnant women or women who are planning to become pregnant.

## **POTENTIAL ADVERSE EVENTS / TECHNICAL COMPLICATIONS**

Possible adverse events include, but are not limited to, the following:

- Access site thrombosis
- Air embolism
- Atrial fibrillation
- Blood flow restriction
- Caval thrombosis or caval occlusion
- Damage to the inferior vena cava
- Death
- Deep vein thrombosis
- Embolization of the device resulting in cardiac arrhythmia or compromise of cardiac valve function
- Extravasation of contrast material at time of vena cavogram causing tissue damage
- Filter fracture
- Filter migration
- Filter tilting
- Hematoma or nerve injury at access or conversion site
- Hemorrhage
- Infection, including access site infection
- Perforation of the inferior vena cava and/or adjacent organs or vertebral bodies
- Pulmonary embolism
- Shortness of breath
- Venous insufficiency

## **IMPORTANT INFORMATION, GENERAL WARNINGS AND PRECAUTIONS**

Important information regarding use of the VenaTech® Convertible™ Vena Cava Filter System includes the following:

### **IMPORTANT INFORMATION**

#### **Instructions for Use – Jugular Approach**

These Instructions for Use must be used when implantation of the VenaTech Convertible Vena Cava Filter using the **jugular approach** has been chosen. The Instructions for Use provided in the package for the **femoral approach** must be set aside at the start of the jugular implantation procedure.

#### **Target Population**

This filter is designed to provide effective protection against pulmonary embolism where the patient's risk is temporary or permanent. In cases where the physician determines that it is clinically necessary to discontinue filtration, the intraluminal filter elements of the VenaTech Convertible Vena Cava Filter can be converted to an open configuration.

#### **Access Site**

The VenaTech Convertible Vena Cava Filter is designed to be implanted using the **right internal jugular vein**. Other access routes (including the left internal jugular vein or the right or left external jugular veins) should not be attempted.

**Warning:** Do not use an access site previously used for the implantation of a central venous catheter. The implantation of a vena cava filter using an existing puncture site can result in an incomplete filter deployment due to the presence of tissue encapsulating the filter. This could result in possible filter migration and/or insufficient protection from pulmonary embolism.

**Warning:** Extreme care must be used for approaches with tortuous venous anatomy. Such venous anatomy can cause kinking of the introducer sheath and make filter insertion difficult or impossible. In extreme conditions, attempting to force the filter through a kinked sheath could result in a sheath perforation. When unusual resistance is encountered, stop advancing the filter and remove the filter and introducer sheath assembly.

## Diameter of Vena Cava

Do not deploy the filter unless the IVC diameter has been properly determined to avoid filter migration. The diameter of the vena cava, as evaluated during cavography, must be less than or equal to ( $\leq$ ) 28 mm from an anteroposterior (AP) projection when corrected for magnification. In cases where measurement is accomplished using Computed Tomography (CT), the section with the largest diameter should be used.

**Warning:** Do not attempt to increase the diameter of the filter; the integrity of the filter will be compromised.

## Filter Placement In Case Of Pregnancy

The VenaTech® Convertible™ Vena Cava Filter is contraindicated for placement in pregnant women or women who are planning to become pregnant; female patients should use appropriate birth control measures.

## Choosing the Site for Deployment

Whenever possible, place the filter just below the level of the lowest renal vein. Place the filter so the top of the cone is 0.5 cm below the level of the lowest renal vein.

## Procedure for Conversion

In cases where the physician determines that it is clinically indicated to discontinue filtration, the intra-luminal filter elements of the VenaTech Convertible Vena Cava Filter can be converted to an open configuration. The time period to conversion for clinical study subjects ranged from 15 days to 391 days.

These Instructions for Use must be kept for guidance during the conversion procedure. A current electronic version of the Instructions for Use may also be found at [www.bisusa.org](http://www.bisusa.org).

**Warning:** Conversion of the VenaTech Convertible Vena Cava Filter should not be performed if thrombus has been captured in the filter.

**Warning:** It is necessary to perform a vena cavogram in order to confirm the patency of the IVC and to ensure that there is no further risk of pulmonary embolism before converting the filter.

**Warning:** Conversion should not be performed if the filter is tilted or if the filter has been deployed asymmetrically. Attempting to convert the filter in this instance could result in damage to the filter or to the inferior vena cava (IVC).

## GENERAL WARNINGS

- The VenaTech® Convertible™ Vena Cava Filter system should only be used by physicians who have experience with the technical and clinical requirements of vena cava filter implantation.
- The VenaTech Convertible Vena Cava Filter is not a retrievable filter. Attempting to retrieve the filter will result in damage to the filter or to the inferior vena cava (IVC).
- **STERILE - DO NOT RESTERILIZE - FOR SINGLE-USE ONLY.**  
The VenaTech Convertible Vena Cava Filter System is provided STERILE and intended for single-use only; do not resterilize, reuse, or clean. The device and its accessories are not reusable nor designed to be re-used or cleaned. Any re-use or cleaning may compromise the performance and safety of the device.
- Do not deploy the filter unless the IVC diameter has been properly determined to avoid filter migration.
- Never advance the guidewire or introducer catheter or deploy the filter without fluoroscopic guidance to avoid any introduction mistake (e.g., level of deployment, deployment in collateral veins).
- Do not directly manipulate the implanted filter. Migration of vena cava filters has been reported following direct manipulation of the IVC or its compression during abdominal surgery and, especially, retroperitoneal surgery.
- Migration of the filter could occur during specific activities. Patient should avoid strenuous physical exercises, which increase the intra-abdominal pressure (e.g., Valsalva maneuver).
- After use, the VenaTech Convertible Vena Cava Filter System and accessories may be a potential biohazard. Handle and dispose of these materials in accordance with current medical practices and applicable laws and/or regulations.
- The VenaTech Convertible Vena Cava Filter system should not be used on patients with known allergies to the following materials: Cobalt Chromium metal alloy (CoCr20Ni16Mo7), Polymethylmetacrylate, Polyethylene, Polypropylene, Polycarbonate, Fluorinated ethylene propylene, stainless steel.

## GENERAL PRECAUTIONS

- Secondary catheterization: Following filter placement, no central venous catheterization should be attempted without fluoroscopic guidance to ensure the location of the deployed filter is known thereby reducing potential for entrapment of other devices in the filter and subsequent sequelae.
- Verify that the sterile package is undamaged before opening. Do not use if the package is damaged.
- Confirm the date of sterile expiration on the labeling. Never implant a VenaTech Convertible Vena Cava Filter that has exceeded its expiration date. To maintain sterility, all normal operating room sterile procedures should be followed.

## MRI SAFETY INFORMATION

### MRI Safety Information

Non-clinical testing has demonstrated the VenaTech® Convertible™ Vena Cava Filter is MR Conditional. A patient with this device can be safely scanned in an MR system meeting the following conditions:

- Static magnetic field of 1.5T or 3T
- Maximum spatial field gradient of 1,160 gauss/cm (11.6 T/m)
- Maximum MR system reported, whole body averaged specific absorption rate (SAR) of 2 W/kg (Normal Operating Mode)

**Note:** The MRI information is related to the implanted filter only. The implantation accessories were not tested.

Under the scan conditions defined above, the VenaTech Convertible is expected to produce a maximum temperature rise of less than 3.3° C after 15 minutes of continuous scanning.

In non-clinical testing, the image artifact caused by the device extends approximately 11 mm from the VenaTech Convertible when imaged with a gradient echo pulse sequence and a 3T MRI system.

## **CLINICAL STUDY SUMMARY**

Study Design: A U.S. multi-center, prospective, single-arm, clinical study was conducted to assess the safety, performance, and convertibility of the VenaTech® Convertible™ Vena Cava Filter.

The primary objective of the study was to demonstrate that the rate of technical success is no lower than the rate (retrieval success) reported in the literature control. Technical success was defined as filter conversion without the loss of filter head components in the vasculature or incomplete opening of filtering legs. The objective performance criterion (OPC) for technical success was 80%. During the analysis, the sponsor did not count any filters as a 'technically' successful conversion when the operator was unable to snare the filter hook during an attempted conversion.

The secondary objective was to estimate the 6-month major device-related adverse event rate in subjects with a converted filter. Major adverse events are defined as:

- Pre-conversion: Documented symptomatic pulmonary embolism, symptomatic caval thrombosis, clinically significant filter migration, or IVC perforation.
- Post-conversion: Symptomatic caval thrombosis, clinically significant filter migration or IVC perforation.

The ancillary measure was to determine the device deployment success rate in subjects where there was an attempt to place a filter. The device was considered successfully deployed when it was advanced and implanted at its intended location in the IVC without:

- A delivery system failure
- Embolization of filter
- Extravascular perforation of IVC by guidewire or filter components
- Incomplete opening of stabilizing legs
- Inadequate distribution of filtering legs (i.e., filtering mechanism)
- Misalignment with axis of IVC (e.g., significant tilt)
- Significant filter migration

A total of 149 subjects underwent filter placement (80 male, 69 female; mean age 62.7 years; age range, 19.8 – 90.3 years) at 11 sites across the United States. All subjects were followed for 6-months post-implant including a clinical exam and KUB and Doppler or KUB and CT imaging. The standard of care baseline, implant procedure, and protocol-required 6-month radiographic imaging was reviewed by a central imaging core laboratory. In addition to the 6-month visit, converted subjects were also followed at 30-day and 3-months post filter conversion to assess for adverse events. All reported adverse events were reviewed and adjudicated by an independent physician Clinical Events Committee.

Indications for filter implant were recorded on the case report form by the investigators as presented in Table 1.

Table 1: Indication for Filter Placement	
Condition/Risk Factor <sup>1</sup>	All Subjects (n=149)
Pulmonary thromboembolism when anticoagulants are contraindicated (%)	20.8% (31/149)
Failure of anticoagulant therapy in thromboembolic diseases (%)	12.1% (18/149)
Emergency treatment following massive pulmonary embolism where anticipated benefits of conventional therapy are reduced (%)	7.4% (11/149)
Chronic, recurrent pulmonary embolism where anticoagulant therapy has failed or is contraindicated (%)	0.7% (1/149)
Other (%)	59.1% (88/149)
<sup>1</sup> <b>Note:</b> Sites had the option to choose more than one indication for IVC filter placement, therefore, the sum of the numerators will not always add up to the denominator.	

The majority of responses were categorized by the investigators as 'other', suggesting the investigators did not think the actual indication for filter placement was adequately described in one or more of the pre-defined categories listed on the case report form. All 'other' responses were further categorized into groupings discussed with the CEC, as follows:

- Acute pulmonary embolism (1 subject)
- Current DVT (2 subjects)
- Current DVT with contraindication to anticoagulation (24 subjects)
- Pre/post surgery prophylaxis (16 subjects)
- Pre/post surgery prophylaxis with remote history of PE/DVT (31 subjects)
- Trauma (14 subjects)

Subject accountability is presented in Table 2.

Table 2: Accountability					
Time Point	Eligible for Visit (N)	Visit Completed (N, %)	Death (N)	Lost to Follow-Up (N)	Withdrawal (N)
Baseline / Implant	149	149, 100%	0	0	0
Permanent Filtration 6-Month	40	39, 97.5%	11	1	3
Converted	96	96, 100%	0	0	0
Converted 30-Day	92	91, 98.9%	1	1	0
Converted 3-Month	91	89, 97.8%	1	0	1
Converted 6-Month	90	89, 98.9%	0	0	0

Results: The protocol ancillary measure was to determine the device deployment success rate in subjects where there was an attempt to place the filter. All filter implant attempts were successful (100% 149/149) with no reports of filter malposition. Of the 149 enrolled subjects, 121 were assessed for filter conversion eligibility. Filters remained unconverted for a variety of reasons, similar to the reasons retrievable filters are not retrieved, including subject withdrawals due to death, physician or subject preference, continued risk of pulmonary embolism, withdrawal of consent, or lost to follow-up. The primary objective was met demonstrating the technical success of filter conversion exceeded the OPC of 80%. There have been 96 filter conversion attempts reported, 89 of which have no associated technical complications affecting the primary endpoint, resulting in a 92.7% conversion success rate (lower 97.5% confidence limit 85.6%). There were no reports of loss of filter head components in the vasculature, four (4) reports of incomplete opening of the filtering legs immediately post conversion after use of ancillary tools, and three (3) reports of inability to convert the filter due to the user being unable to snare the filter hook. Details of the technical complications associated with conversion procedures are as follows.

#### **Converted with incomplete opening of the filter legs (4):**

- In subject 001-003, the investigator was unable to disengage the filter cap easily (81 days post-implant). The investigator used extra maneuvers to disengage the cap, which involved a second access from the subject's groin. The investigator used a cobra catheter to push the filter cap from below while pulling the cap from above at the jugular access. Once the cap was disengaged, accessories were used to convert the filtering legs to an open configuration (balloon inflation and reverse curve catheter). The "majority" of the legs deployed into the expected "stent" (open) configuration. (Note: the core lab confirmed the investigator's assessment and said the filter legs are tethered together but the core lab could not tell if the legs were apposed to the IVC because of the cavogram view.) No thrombus or filling defect was identified in the filter or elsewhere in the IVC. No adverse events were reported associated with this conversion procedure.
- In subject 004-004, the investigator converted the filter (158 days post implant). According to the investigator, "the head was easy to snare and release but was difficult to disengage from the filter. The struts did not open at all. A cobra catheter, tip deflecting wire, and balloon were used to partially open the struts." The core lab agreed with the investigator's assessment confirming the legs did not open completely. No adverse events were reported associated with this conversion procedure.
- In subject 005-001, the investigator converted the filter 169 days post implant. The filter hook was successfully snared and the filter head components successfully retrieved, however all filter legs did not completely open to a stent like position even after use of an SOS catheter. Per the core lab's review of the 6-month imaging (CT), the technical complication is resolved and the legs are fully open (apposed to IVC wall). No adverse events were reported associated with this conversion procedure.

- In subject 007-019, the investigator converted the filter 207 days post implant. The filter hook was successfully captured and the filter head components successfully retrieved, however all filter legs did not completely open into a stent-like position even after use of a variety of ancillary tools including a compliant balloon, RIM catheter, and deflecting tip wire, which was used to attempt to try to free three of the legs. The remaining legs appeared to be embedded within the wall of the IVC. The core lab reviewed the subject's 6-month post-conversion follow-up visit imaging and noted "the filter is nearly at full IVC wall apposition on latest images". No adverse events were reported associated with this conversion procedure.


### Inability to convert the filter (3):

- In subject 001-016, the investigator attempted to convert the filter (180 days post implant) and was unable to snare the hook. The filter remains unconverted. Neither the site nor the core lab reported significant filter tilt or misalignment. The absolute cause of the technical complication (inability to snare and subsequently convert filter) is unknown but may be due to hook endothelialization. A conversion-related adverse event of 'back pain' was reported.
- In subject 007-011, the investigator attempted to convert the filter (159 days post implant) and was unable to snare the filter hook due to investigator reported hook endothelialization (however, the core lab reported there was no evidence of endothelialization based on the images received). Neither the investigator nor the core lab reported filter tilt. The filter remains unconverted. No adverse events were reported associated with this conversion attempt.
- In subject 009-016, the investigator attempted to convert the filter (212 days post implant) and was unable to snare the filter hook; the filter remains unconverted. Neither the site nor the core lab reported significant filter tilt or misalignment. No adverse events were reported associated with this conversion attempt.

### Five (5) conversion procedures were classified by investigators as 'difficult', as follows:

- Subject 001-003 – Inability to easily disengage the filter cap.
- Subject 001-019 – Per investigator, due to angulation (7 degrees) of filter, the hook had endothelialized and had to be moved from the caval wall with a balloon so the filter head could be removed.
- Subject 004-004 – Per investigator, the head was easy to snare and release but it was difficult to disengage from the filter. The struts did not spontaneously open. A cobra catheter, tip deflecting wire, and balloon were used to partially open the struts of the filter.
- Subject 006-003 – Per investigator, fibrin sheath was present on the filter requiring multiple maneuvers to open filter legs.
- Subject 007-019 – Inability to secure full wall apposition of filtering legs despite use of multiple accessories.

The mean days to filter conversion was 130.7 days (range 15 to 391 days). The average reported conversion procedure time was 30.7 minutes (the conversion procedure duration was defined as the time the introducer is placed to the time all ancillary accessories were removed). In 82.3% of conversion procedures, accessory tools were used to assist including guiding catheters, guidewires and angioplasty balloons (as allowed per the instructions for use). Figure 1 presents days to conversion attempt.

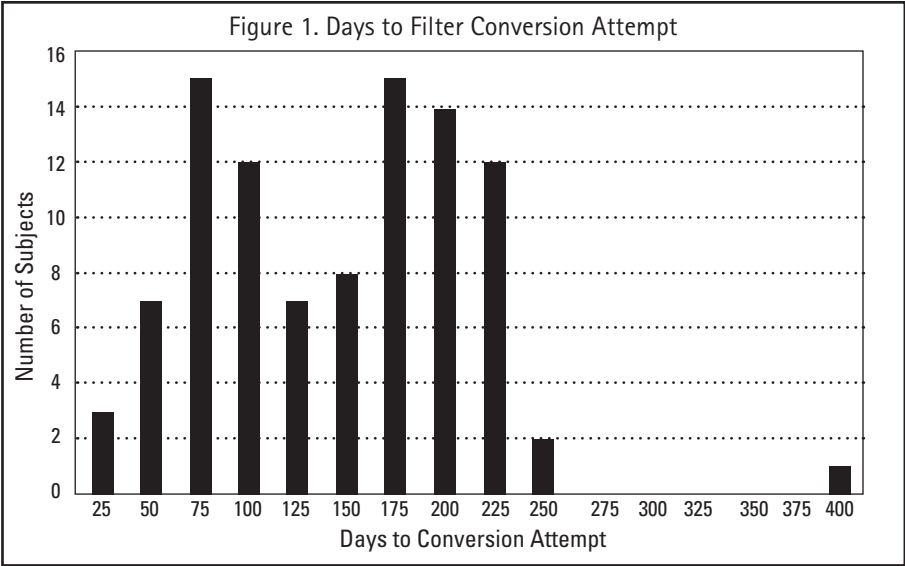**Adverse Events & Technical Complications:** The literature supported expected rate for the secondary endpoint of major device-related adverse event rate is 3 – 6 %. There were no major device-related adverse events reported for converted subjects through 6-months. The estimated rate of device related major adverse events is therefore 0.0% with 95% binomial exact confidence limits of 0.0% to 4.1%.

Table 3 and Table 4 present reported technical complications and relevant adverse events. There have been no reports of device or procedure-related major adverse events post-conversion (defined as symptomatic caval thrombosis or caval occlusion, perforation of the IVC and/or adjacent organs or vertebral bodies, pulmonary embolism, or filter migration). There have been no reports of spontaneous filter conversion or loss of the filter head during the conversion procedure.

In summary, the results of this study demonstrated that the placement and conversion of the VenaTech® Convertible™ Vena Cava Filter can be performed safely with relatively high rates of procedural and clinical success. For patients who are no longer at risk for thromboembolism, the VenaTech Convertible filter can be implanted for several months and then safely converted. These data demonstrates the safety and effectiveness of the placement and conversion of the VenaTech Convertible filter system in a clinically relevant patient population when used according to the instructions for use.

**Table 3: All Reported Technical Complications**  
(Shaded rows represent no reported event)

Technical Complications	All (n=149)
Deployed at unintended position	0.0% (0)
Embolization of filter or filter components	0.7% (1) <sup>1</sup>
Filter fracture	0.0% (0)
Filter implanted upside down	0.0% (0)
Filter migration, significant (> 20 mm)	0.7% (1) <sup>1</sup>
Inadequate distribution of filtering legs	0.0% (0)
Incomplete opening of filtering legs during conversion procedure	2.7% (4)
Incomplete opening of stabilizing legs during deployment	0.0% (0)
Misalignment with axis of IVC / tilt (> 15°)	0.7% (1)
Spontaneous conversion of filter	0.0% (0)
Other	4.0% (6) <sup>2</sup>

<sup>1</sup> One un-converted subject experienced an unanticipated adverse device effect (UADE), significant filter migration resulting in atrial perforation. The DSMB believes this event was related to the "sail effect" of a large clot burden reaching the filter and carrying the filter with it to the heart, based on the large clot in the filter, and the associated significant PE. This event resulted in repeat surgical explant procedures to retrieve filter remnants. Of the total eight related serious adverse events reported in the study, seven occurred in this subject. Complications include pulmonary embolism, right atrial thrombus, right ventricle thrombus, right ventricle perforation, compromise of cardiac valve function due to filter embolization, mitral valve regurgitation and valve repair.

<sup>2</sup> Other reported technical complications:

- Unable to convert filter (n=3).
- Post conversion cavogram demonstrated a small filling defect along the lateral wall of the IVC at the superior margin of the device, likely representing a small dissection flap – not flow limiting per investigator. A balloon was reinserted and inflated (n=1).
- One left sided filter strut tilted toward center of IVC; identified by core laboratory on 6-month imaging (n=1).
- One report of IVC penetration < 3 mm (n=1). IVC penetration (a chronic process) defined as post-implant penetration of the IVC wall by one or more filter components that extends > 3 mm outside the wall of the IVC as documented by cross-sectional imaging when there is a clinical or radiographic suggestion that penetration may have occurred.

**Note:** This study distinguished between IVC penetration and IVC perforation. IVC extravascular perforation (an acute process) is defined as during placement of the filter, a guidewire or filter component pierces and extends outside the IVC wall by more than 3 mm as documented by radiographic imaging. There were no reports of IVC perforation (see Table 4).

**Table 4: Relevant Adverse Events**  
(Shaded rows represent no reported event)

Anticipated Events	Overall		Device Related		Implant Procedure Related		Conversion Procedure Related	
	Events N	% of Subjects (n/N)	Events N	% of Subjects (n/N)	Events N	% of Subjects (n/N)	Events N	% of Subjects (n/N)
Access site thrombosis	1	0.7% (1/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Air embolism	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Cardiac arrhythmia due to filter embolization	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Caval thrombosis or caval occlusion	5	3.4% (5/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Compromise of cardiac valve function due to filter embolization	1	0.7% (1/149)	1 <sup>1</sup>	0.7% (1/149)	0	0.0% (0/149)	0	0.0% (0/96)
Damage to inferior vena cava	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Death	13	8.7% (13/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Deep vein thrombosis (DVT)	23	10.1% (15/149)	0	0.0% (0/149)	1	0.7% (1/149)	0	0.0% (0/96)
Hematoma at access site or conversion site	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Infection, including access site infection	4	2.0% (3/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Nerve injury at access site or conversion site	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
IVC perforation and/or perforation of adjacent organs or vertebral bodies > 3 mm <sup>2</sup>	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Pulmonary embolism	1	0.7% (1/149)	1 <sup>1</sup>	0.7% (1/149)	0	0.0% (0/149)	0	0.0% (0/96)
Tissue damage caused by extravasation of contrast material at time of venacavogram	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Venous insufficiency	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/149)	0	0.0% (0/96)
Other			6 <sup>1,2</sup>	1.3% (2/149)	2 <sup>3</sup>	1.3% (2/149)	3 <sup>4</sup>	3.1% (3/96)

<sup>1</sup> Five of six device-related 'other' adverse events occurred in the subject described above in Table 1, footnote 1.

<sup>2</sup> Thrombus above the filter (n=1).

<sup>3</sup> Loss of vascular access during filter implant (n=1) and minor bruising at insertion site without hematoma (n=1).

<sup>4</sup> Back pain (n=1), epidermoid cyst at conversion access site (n=1), and soreness at access site post-conversion (n=1).

## **IMPLANTATION PROCEDURE**

Prior to implanting the VenaTech® Convertible™ Vena Cava Filter, a pre-placement cavogram is required to:

- Confirm the patency of the vena cava.
- Evaluate the anatomy of the vena cava and identify any venous anomalies.
- Mark the level of the renal veins.
- Determine the highest level of any thrombus which may be present.
- Decide on the desired level for filter deployment and to mark the intended position with respect to the vertebral bodies.
- Ensure the diameter of the vena cava (AP projection) at the planned filter deployment site is less than or equal to ( $\leq$ ) 28 mm in diameter.

### **A) Preparation of the placement accessories**

Immerse the pre-loaded filter cartridge in heparinized saline.

**Note:** Immersion in heparinized saline reduces the possibility of thrombus formation around the filter during its movement through the introducer sheath.

### **B) Puncturing the vein**

1. Prepare the patient for access using the jugular approach.
2. Echographic puncture the vein using the Seldinger technique (or using cut-down). Confirm venous reflux in the syringe.
3. Advance a 0.035" guidewire and then remove the Seldinger needle.  
(**Note:** This guidewire is not provided in the VenaTech Convertible Vena Cava Filter System package.)

### **C) Pre-placement control cavogram and marking the filter implant site.**

1. Make a small incision using a scalpel to enlarge the access site and facilitate the passage of the introducer sheath. Insert a 9F catheter introducer sheath over the guidewire. Remove 'J' guidewire. Insert angiography catheter (or catheter of your choice) and perform a retrograde cavogram.

**Warning:** If ilio-caval clot is present, extreme care must be exercised to avoid dislodging the clot.

2. Locate the renal veins and choose the filter deployment site. Mark the position of the top of the filter cone by placing a metallic marker (hemostat) transversely across the patient's abdomen.

**Important Note:** Choosing the site for deployment: Whenever possible, place the filter just below the level of the lowest renal vein. Place the filter so that the top of the cone is 0.5 cm below the level of the lowest renal vein.

#### D) Insertion of a coated 'J' guidewire


Remove the angiographic catheter and insert a coated 'J' guidewire.

**Warning:** To prevent the introducer from inadvertently passing outside of the IVC (e.g., into a tributary vessel or the right ventricle) verify that the 'J' portion of the guidewire reaches the iliac bifurcation.

#### E) Placement of the filter introducer sheath

1. Remove the 9F catheter introducer sheath.
2. Insert and advance the filter dilator/sheath combination over the 'J' guidewire.
3. Continue to advance the filter dilator/sheath combination until the introducer sheath marker is slightly more than 6 cm below the filter deployment site marker (hemostat) placed on the abdomen.

**Note:** The marker on the introducer sheath notes the level of the bottom of the filter when it is within the sheath just prior to deployment.

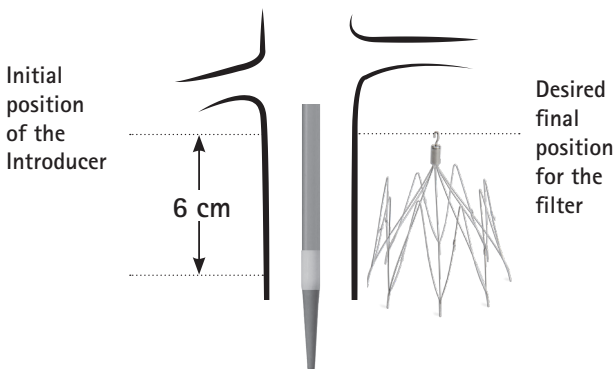This maneuver pre-positions the filter slightly lower than the level chosen for its deployment. The introducer system is drawn upward slightly prior to deployment to adjust filter level to its final desired position.

**Note:** The final position adjustment should occur just prior to filter deployment.


**Warning:** Hold the guidewire firmly while advancing the introducer to prevent its movement or kinking.

## F) Inserting the filter in the sheath and advancing it into position

1. Disconnect the luer-lock securing the dilator to the sheath.
2. Remove the 'J' guidewire and the dilator together taking care not to move the sheath.
3. Flush the sheath with 10 ml of heparinized saline.
4. Connect the filter cartridge securely to the sheath (see image):

5. Using the pusher, advance the filter through the sheath under fluoroscopy until the proximal marker on the pusher reaches the top of the filter cartridge.

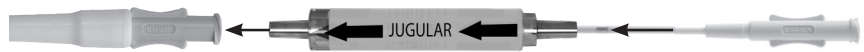
**Warning:** The cartridge must remain connected to the sheath while using the pusher in order to avoid premature filter ejection.

6. The filter is now positioned at the distal end of the sheath. The system is now ready for final deployment positioning.

## G) Filter deployment

1. Under fluoroscopy, adjust the sheath position so that the top of the filter cone (visible thru the sheath under fluoroscopy) is at the level of the radiopaque marker (hemostat) placed on the abdomen of the patient at the beginning of the procedure. Ensure that the top of the filter cone is 0.5 cm below the level of the lowest renal vein. This is the final position adjustment for the filter prior to deployment.
2. When the filter is properly positioned, avoid any further movement of the introducer system.
3. Confirm using fluoroscopy that the sheath has not moved and that the filter is properly positioned at the desired level for placement.

4. Holding the pusher in one hand, take the sheath in the other hand and under fluoroscopy, draw the sheath towards yourself. This movement deploys the filter.

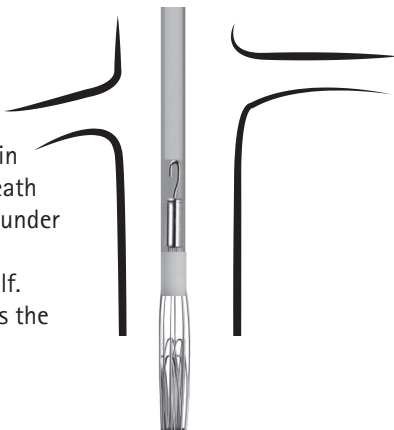
**REMINDER:**

To ensure proper filter positioning, it is extremely important to not move the pusher:

- After the filter is positioned for deployment.
- During filter deployment.

5. The filter is now implanted

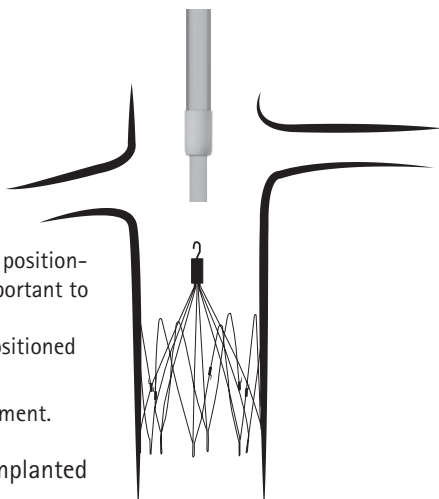**Note:** To perform a control cavogram at this point in the procedure:

- Withdraw the pusher/sheath/cartridge assembly several centimeters.
- Remove the pusher and the cartridge together while keeping the sheath in place.
- Perform the cavogram by manually injecting contrast material through the sheath, then remove the sheath.

6. Remove the sheath assembly.

## H) Final documentation

Take a final film KUB (kidneys, urethra, bladder), which documents the location and deployment of the filter.

## **CONVERSION PROCEDURE**

In cases where the physician determines that after filter implantation it is clinically indicated to discontinue filtration, the VenaTech® Convertible™ Vena Cava Filter can be converted to an open configuration to no longer function as a filter. Filter conversion is achieved by percutaneously unlocking the mechanism that secures the filter head to the filtering legs and removing the filter head. The time period to conversion for clinical study subjects ranged from 15 days to 391 days.

### **REQUIRED MATERIALS FOR THE CONVERSION PROCEDURE:**

These materials are not provided in the VenaTech Convertible Vena Cava Filter System package.

Physician must have the following items available to facilitate conversion and ensure filter legs are fully apposed to vessel wall:

- 'J' Guidewire (0.035 in / 0.89 mm, 150 cm length)
- Catheter introducer sheath (9F I.D.) (to facilitate removal of the filter head)
- Amplatz GooseNeck™ snare (sized matched to the diameter of the vena cava) and catheter
  - 15 mm loop and 6F catheter (I.D. 1.6 mm), length ≥ 102 cm
  - 20 mm loop and 6F catheter (I.D. 1.6 mm), length ≥ 102 cm
  - 25 mm loop and 6F catheter (I.D. 1.6 mm), length ≥ 102 cm
- Accessory device such as a diagnostic catheter, guidewire, or angioplasty balloon

#### **A) Perform vena cavogram**

Perform a vena cavogram as described in the Implantation Procedure in this IFU to confirm the absence of thrombus in or below the filter (i.e., peripheral vasculature).

**Warning:** Conversion of the VenaTech Convertible Vena Cava Filter should only be performed using an Amplatz GooseNeck™ snare. Use of another snare may result in complications, including embolization of the filter head and/or snare fracture.

**Warning:** Conversion of the VenaTech Convertible Vena Cava Filter should not be performed if thrombus has been captured in the filter.

**Warning:** It is necessary to perform a vena cavogram in order to confirm the patency of the IVC and to ensure that there is no further risk of pulmonary embolism before converting the filter.

**Warning:** Lack of traction on the snare may lead to the loss of the filter head in the circulatory system.


**Warning:** Optimal conversion is facilitated by having the external parts of the snare catheter and the snare in the same plane as the patient which provides minimal torque loss and the most predictable snare response at the distal end.

## **B) Placement of the intravascular snare**

1. Select the approach which provides the best access to the filter head for conversion.
2. Confirm that the hook on the cone of the filter is oriented towards the puncture site chosen for filter conversion.
3. Place a 9F catheter introducer sheath using the Seldinger technique.
4. Position the 'J' guidewire (150 cm length) in the inferior vena cava.
5. Advance the snare catheter over the 'J' guidewire to the top of the hook on the filter cone.
6. Remove the 'J' guidewire.
7. Introduce the snare through the snare catheter to the top of the hook on the filter cone.

## **C) Capture the hook**

1. Capture the hook on the cone of the filter with the snare loop.
2. Exert tension on the snare and straighten snare completely.
3. Advance the snare catheter over the snare to ensure that the snare catheter is 1 cm from the hook while maintaining tension on the snare.
4. Ensure that the snare catheter, the snare, and the hook are in the same plane.

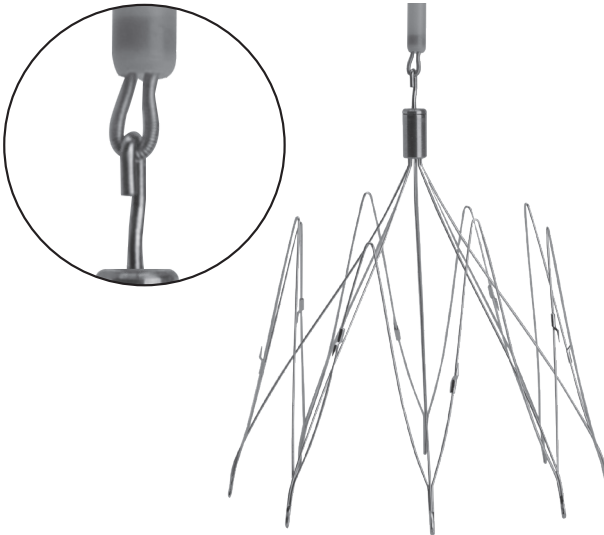
## D) Verification of snare and catheter position

1. Maintain tension on the snare and advance the snare catheter downwards until it covers the hook.

### HEAD LOCKED:

The snare catheter radiopaque marker is visible between the top of the hook and the head of the filter

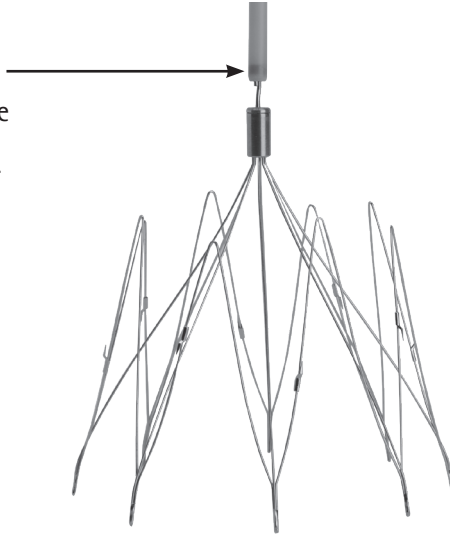2. Confirm that the profile of the hook and filter head conforms to the above corresponding image.

## E) Unlocking the filter head

1. Keep the proximal part of the snare catheter and snare straight while under tension.
2. Under fluoroscopy, maintain the snare catheter position as referenced in the previous image. Pin the catheter in place and pull the snare proximally until the filter head unlocks.

**Note:** If the filter head does not unlock immediately start the procedure again from Step D (Verification of snare and catheter position) above.

3. Confirm that the profile of the hook and the head corresponds to the following image.

### HEAD UNLOCKED:

The marker inside the snare head is in its upward position. The distance between the snare catheter radiopaque marker and the top of the hook has increased.

## F) Retrieval of the head

1. Maintain traction on the snare to ensure that the hook remains in the snare catheter at all times.

**Warning:** Lack of traction on the snare may lead to the loss of the filter head in the circulatory system.

**Warning:** In case of difficulty retrieving filter head, do not apply excessive force to the filter.

2. Remove the snare catheter, snare, and filter head as a single unit.

## **FINAL DOCUMENTATION**

Take a final cavogram and a final film (KUB) to verify that the filter legs have completely opened and are in apposition with the vena cava wall (see image below). (Incomplete opening of the filter legs is defined as less than complete apposition of the legs to the vena cava wall as evidenced by cavography.)

### **'Converted' Filter (filter head removed)**

**Note:** In the event that the filter does not fully open immediately after the filter head is removed, an accessory device may be used to complete the filter conversion. See below for instructions regarding filter conversion with accessories.

## **CONVERSION WITH ACCESSORIES**

After filter head removal, use of an accessory to aid in the filter conversion is probable due to the fibrin strands or cellular growth that collect at the top of the filtration cone. It is recommended that the physician be prepared to utilize an accessory device such as a diagnostic catheter, guidewire, or angioplasty balloon to assist in completely opening constrained filter legs. Insert and advance the accessory device into the introducer and gently coax the filtering leg(s) until they are free from constraining tissue and opened completely. Remove the accessory device.

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