



SECURELY AND CONFIDENTLY
GRASP OBJECTS FROM ANY ANGLE

Multi-Snare®

Dual-Plane Retrieval System


Retrieval and Manipulation of Foreign Objects
from the Vascular System and Hollow Viscera

Multi-Snare®

Dual-Plane Retrieval System

The diameter of the loop on the Multi-Snare® can be easily adjusted to the exact size of the object to be retrieved. In addition to the main loop, the Multi-Snare® also features an additional, lateral retrieval function. The dual plane design enables objects to be grasped quickly, and safely retrieved.

- Low-Profile Introducer
- Dual-plane system increases grasping capability by offering primary and lateral loop catching possibilities
- Variable loop diameter for precise sizing of loop
- Radiopaque loop, platinum coiled with 360° rotation capability
- Flexible nitinol core wire provides high tensile strength
- Low friction with 1:1 torque response

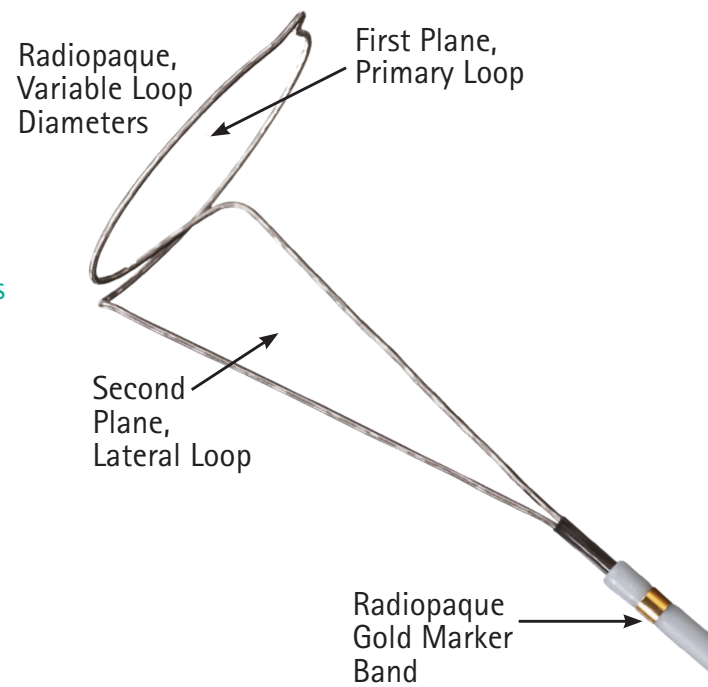

Ordering Information:

Item	REF Number	Re-order Number	Nominal Diameter	Variable Diameter	Snare Length	Recommended Introducer	Introducer Length
Multi-Snare Micro Set	147302	PFM-147302	2mm	2-3mm	175cm	3F	150cm
Multi-Snare Micro Set	147304	PFM-147304	4mm	4-6mm	175cm	3F	150cm
Multi-Snare Set	147305	PFM-147305	5mm	5-8mm	125cm	4F	105cm
Multi-Snare Set	147310	PFM-147310	10mm	10-15mm	125cm	4F	105cm
Multi-Snare Set	147315	PFM-147315	15mm	15-20mm	125cm	5F	105cm
Multi-Snare Set	147320	PFM-147320	20mm	20-30mm	125cm	6F	105cm
Multi-Snare Set	147330	PFM-147330	30mm	30-40mm	125cm	6F	105cm

Multi-Snare® and Multi-Snare® Micro are packaged as a set including Snare, Snare Catheter, Torquer, and Snare Introducer Aid

Indications for Use:


- Retrieval and manipulation of foreign objects from the vascular system and hollow viscera
- Assistance in creating loops where cross-over technique is applied
- Repositioning of indwelling venous catheters
- Assistance in performing venipuncture to obtain access to central vein

For more information or to place an order, contact your B. Braun Interventional Systems Inc. representative or call 1-877-VENA-CAV (836-2228)

Refer to the Multi-Snare® Instructions for Use for relevant warnings, precautions, complications, and contraindications.

LATEX Not made with natural rubber latex

DEHP Not made with DEHP

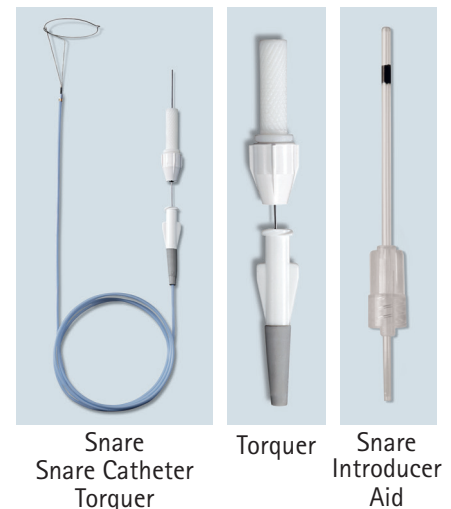

Distributed by:
B. Braun Interventional Systems Inc. | Bethlehem PA | USA
877-836-2228 | www.bisusa.org



Rx only

©2017 B. Braun Interventional Systems Inc.
Multi-Snare is a registered trademark of pfm medical, inc.
CV-331 11/17