

Interventional Systems

**B | BRAUN**

SHARING EXPERTISE



- Hospital ▪ Home
- Nursing Home ▪ Hospice

## ASEPT<sup>®</sup> Pleural and Peritoneal Drainage System

Palliative, Compassionate, and Convenient Patient Care

Needle-Free, Swabable Valve  
Secure Luer Lock Connection  
Active Vacuum Bottle Design

# ASEPT® Pleural and Peritoneal Drainage System

## Catheter Procedure Kit

The ASEPT® Drainage System is a complete insertion kit for pleural or peritoneal catheter implantation. The primary component of the ASEPT® Drainage System is the catheter, which has a self-sealing valve attached to the proximal end that prevents leakage and provides a secure, needle free luer lock connection. The ASEPT® valve is connected to the drainage line on the ASEPT® Evacuated Drainage Bottle or Bag to initiate drainage.

The ASEPT® Drainage Line sets are also available as stand-alone items to allow drainage to compatible drainage collection units.



### TRAY 1:

- Two - ChloraPrep® 10.5mL Tint
- One - Synthetic Suture, 2-0 Straight Needle
- One - Synthetic Suture, 3-0 Curved Needle

### TRAY 2:


- Two - 5cc Syringes
- One - 10cc Syringe
- One - 22 GA. x 1.5" Needle
- One - 25 GA. x 1.0" Needle
- One - 18 GA. Guidewire Introducer Needle
- One - .035" x 70cm J-Tip Guidewire
- One - Safety Scalpel
- One - 16F Tear Away Introducer
- One - 12F Dilator
- One - 8F Dilator
- One - Tunneler
- One - Forceps

## Kit Contents:

- One - ASEPT® Drainage Catheter
- One - CSR Wrap

### TRAY 3:

- One - Fenestrated Drape
- One - 16 GA. x 1.0" Needle
- One - 5-in-1 Drainage Line Adapter
- One - ASEPT® Drainage Line with Protective Cap
- One - Foam Catheter Pad
- Six - 4"x 4" Gauze Pads
- One - 6"x 8" Tegaderm® Self-Adhesive Dressing

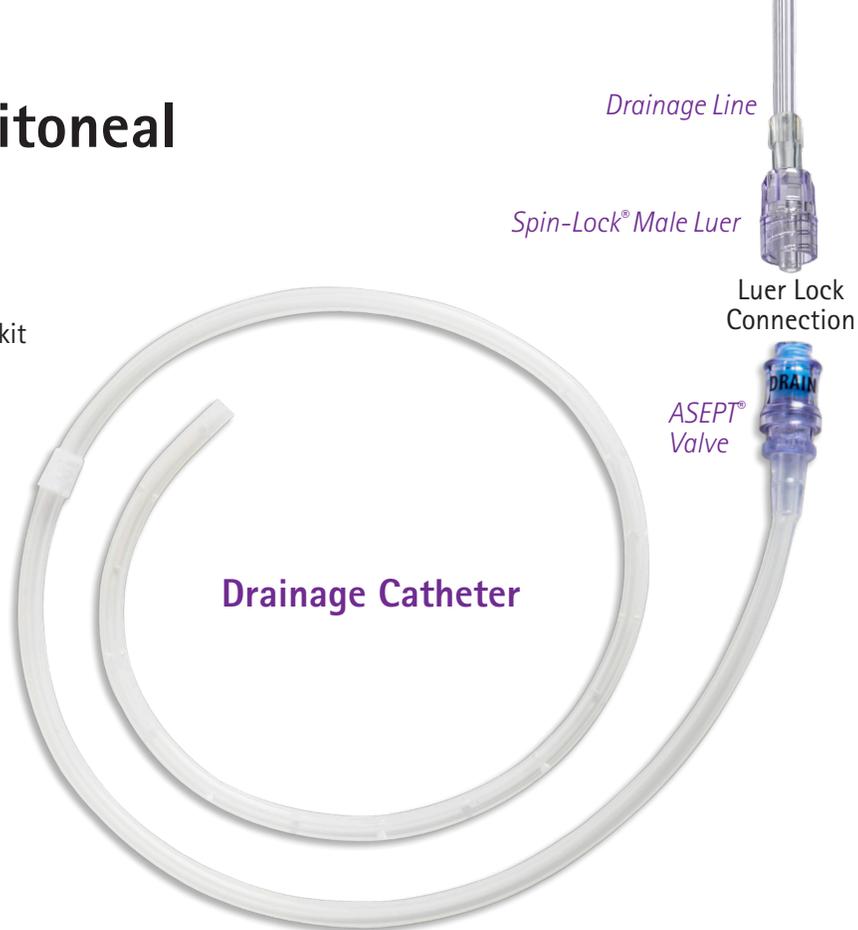

Drainage Line

Spin-Lock® Male Luer

Luer Lock Connection

ASEPT® Valve

Drainage Catheter



# ASEPT® Drainage Kits

600mL and 1,000mL Evacuated Bottles ♦ 2,000mL Bag



## Drainage Procedure Kit

The primary components of the ASEPT® Drainage Kit are the evacuated bottle or drainage bag and the Procedure Pack with all the necessary dressing change items.

The evacuated bottle is available in either a 600mL or 1,000mL configuration. The drainage bag is available in a 2,000mL configuration. The ASEPT® Drainage Kit provides patients with a convenient and compassionate way to perform drainage at home or in any outpatient setting.

### Procedure Pack:

- One - Blue Emergency Slide Clamp
- Three - Alcohol Pads
- One - Foam Catheter Pad
- Six - 4"x 4" Gauze Pads
- Four - Gloves
- One - 6"x 8" Tegaderm® Self-Adhesive Dressing



### Accessories: (Sold Separately)

- ASEPT® Replacement Valve
- ASEPT® Drainage Line Set



## Ordering Information:

ITEM	Reference Number	Case Qty.
<b>ASEPT® Pleural Drainage System</b> (Includes ASEPT® Pleural Drainage Catheter and Insertion Kit)	622289	1
<b>ASEPT® Peritoneal Drainage System</b> (Includes ASEPT® Peritoneal Drainage Catheter and Insertion Kit)	622284	1
		 
<b>ASEPT® 600mL Drainage Kit</b> (Includes ASEPT® 600mL Evacuated Bottle with Drainage Line and Procedure Pack)	622287	10
<b>ASEPT® 1,000mL Drainage Kit</b> (Includes ASEPT® 1,000mL Evacuated Bottle with Drainage Line and Procedure Pack)	622279	10
		
<b>ASEPT® 2,000mL Drainage Bag with Drainage Line</b>	622280	10
<b>ASEPT® 2,000mL Drainage Bag Kit</b> (Includes ASEPT® 2,000mL Drainage Bag with Drainage Line and Procedure Pack)	622281	10
		
<b>ASEPT® Replacement Valve</b>	622288	5
<b>ASEPT® Drainage Line Set</b> (Includes Drainage Line with Protective Cap, Spin-Lock® Luer Connector, 5-in-1 Adapter, and Pinch Clamp)	622286	10
		 

### Indications For Use:

The ASEPT® Pleural Drainage System is intended for long-term, intermittent drainage of symptomatic, recurrent, pleural effusions; including malignant pleural effusions and other pleural effusions that do not respond to treatment of the underlying disease.

The ASEPT® Peritoneal Drainage System is indicated for periodic drainage of recurrent and symptomatic malignant ascites. The catheter is intended for long term access of the peritoneal cavity in order to relieve symptoms such as dyspnea.

The ASEPT® 600mL or 1,000mL Drainage Kit is an accessory to be used to connect to the ASEPT® Pleural or Peritoneal Drainage Catheter to drain fluid from the chest or abdomen.

The ASEPT® 2,000mL Drainage Bag is an accessory to the ASEPT® Peritoneal Drainage System and is used for peritoneal drainage with the ASEPT® Peritoneal Drainage Catheter.

The ASEPT® 2,000mL Drainage Bag Kit is an accessory to the ASEPT® Peritoneal Drainage System and is used for peritoneal drainage with the ASEPT® Peritoneal Drainage Catheter.

The ASEPT® Replacement Valve is used for replacing the valve attached to the ASEPT® Pleural or Peritoneal Drainage Catheter in case of occlusion or malfunction. The ASEPT® Pleural or Peritoneal Drainage Catheter is the main component of the ASEPT® Pleural or Peritoneal Drainage System.

The ASEPT® Drainage Line is an accessory to the ASEPT® Pleural or Peritoneal Drainage System that is used for connecting the ASEPT® Catheter to wall suction or standard fluid collection equipment. The ASEPT® Pleural or Peritoneal Drainage Catheter is the main component of the ASEPT® Pleural or Peritoneal Drainage System.

Refer to the Instructions for Use for complete indications, relevant warnings, precautions, complications, and contraindications.

### Rx only

ASEPT is a registered trademark of pfm medical, inc. Spin-Lock is a registered trademark of B. Braun Medical Inc. Tegaderm is a registered trademark of 3M.

©2018 B. Braun Interventional Systems Inc. Part of the B. Braun Group of Companies.  
CV-1088 Rev. 7/18

Distributed by:

B. Braun Interventional Systems Inc. | 824 Twelfth Avenue | Bethlehem, PA 18018 | USA  
Tel 877-836-2228 | Fax 610-849-1334 | [www.bisusa.org](http://www.bisusa.org)

